


MQ TRI WALKER

mobiquip.com

MQ
MobiQuip

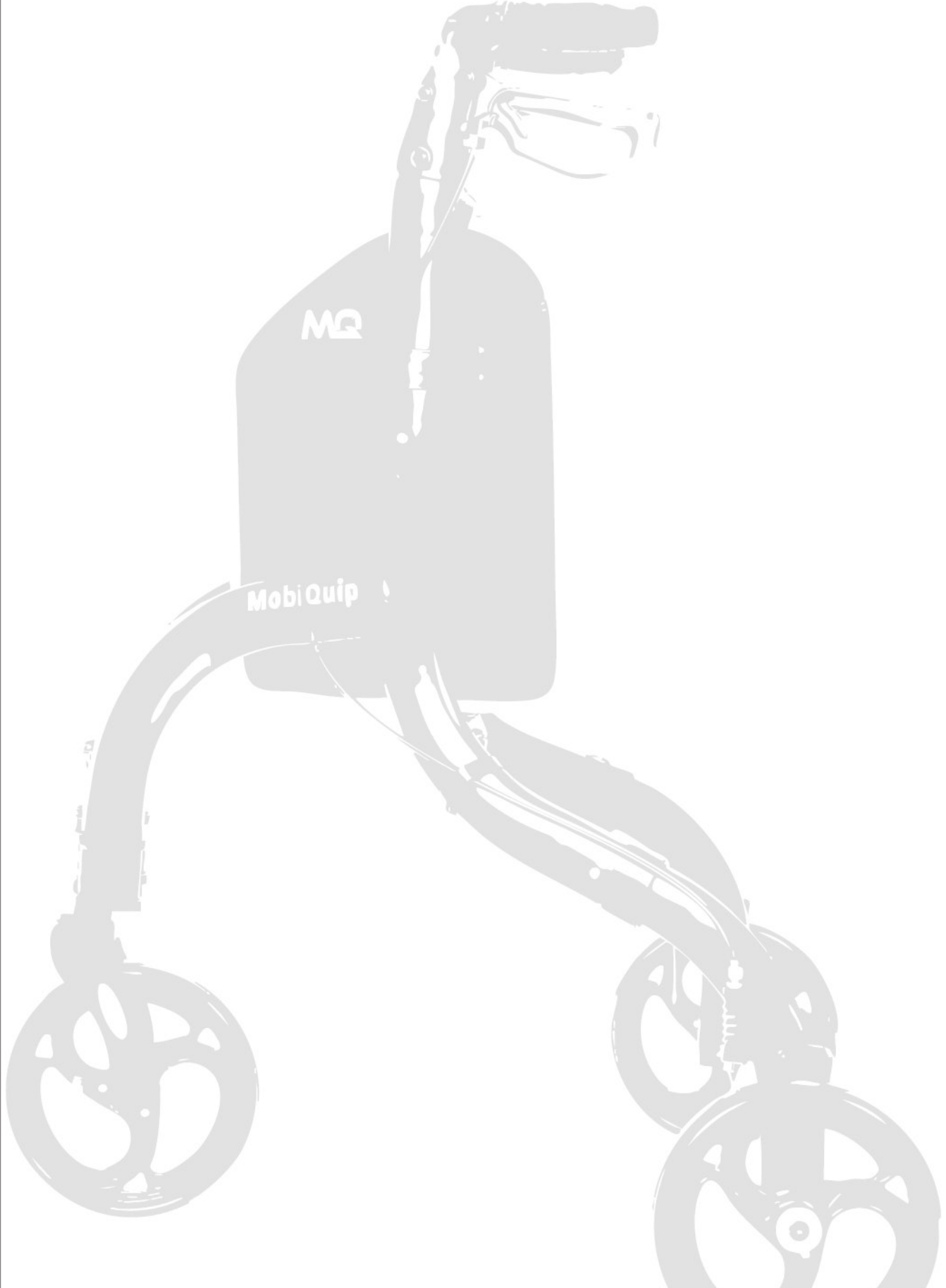

YOU AND/OR YOUR ATTENDANT MUST READ AND FULLY UNDERSTAND THIS MANUAL BEFORE USE.

The MobiQuip Tri-Walker is made of a lightweight aluminium folding frame coated with metallic paint. It is equipped with height-adjustable handles, finished with ergonomic handles and brakes with a parking clamp. The walker has 8" diameter wheels and a zip carry bag.

Features of the MQ Tri-Walker



Application

The walker is designed to support individuals with reduced mobility who require stable assistance while walking. It is particularly beneficial for older adults, helping them maintain balance, walk more comfortably, and reduce the risk of falls. Regular use of the walker can aid in regaining mobility and improving the ability to stand and move independently. It provides enhanced stability, reliable support, and greater confidence while walking.

Assembling, Folding, Unfolding and

Adjustment Assembly elements:



Assembly

- Å Unfold the frame by pulling apart both sections of the crossbar until it locks securely into place (see Fig. 1).
- Å Insert the left rear leg into the left rear socket on the frame. To do this, press the leg locking pin while sliding the leg into place, ensuring the brake wheel is positioned on the outside of the frame. Release the pin so it locks into the hole on the walker frame.
- Å The left leg is connected to the left handle via a brake cable. Insert the left handle into the corresponding hole on the frame, pressing the handle locking pin while sliding it in. Adjust the handle to the desired height, then release the pin so it locks into place.
- Å Ensure the brake cable remains on the outside of the frame and secure it neatly using the Velcro fasteners provided (see Fig. 2).
- Å Repeat the same process for assembling the right rear leg and right handle.




Fitting the Bag

- Insert one of the T-brackets located on the side of the bag into the slot at the top of the main frame (see Fig. 3). Make sure the bag pockets are positioned at the rear of the walker, facing the user.
- Repeat the process on the other side of the bag, securing the second T-bracket into the corresponding slot on the opposite side of the frame.

Brakes

- To stop the walker, pull the brake handle upwards and hold (Fig.4)
- During a stop, lock the parking brake by pushing the brake handle down.
- To release the parking brake, pull the brake handle upwards.

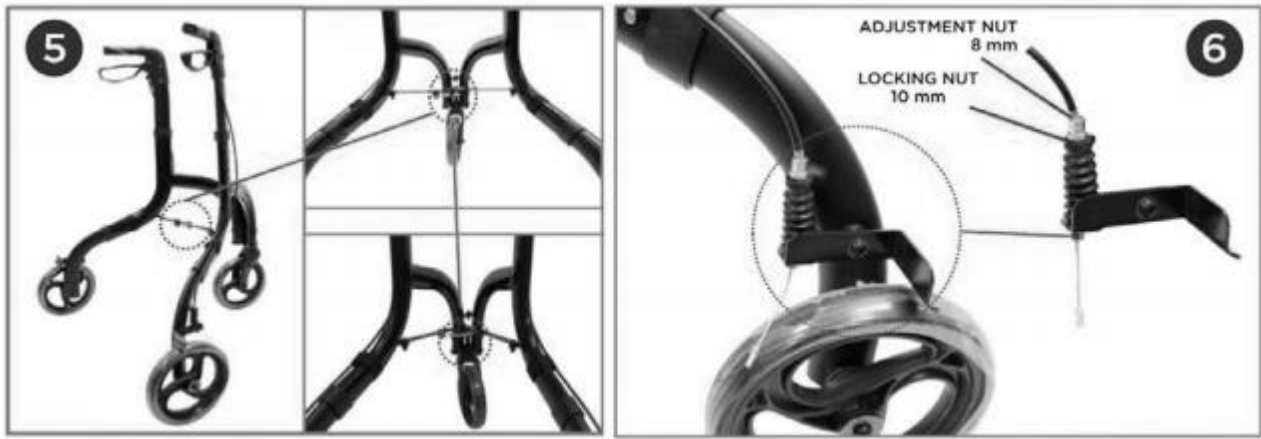


Handles Height Adjustment

- Press both brakes down to set the parking mode and immobilise the rollator.
 - Press and hold the handle locking pin, then adjust the height of the handle.
 - Release the locking pin to lock the adjusted height of the handles. Pull the handle up and down several times to check the correct height locking.
 - Adjust the height of the second handgrip in a similar way.
- Note!** Both handgrips should be set at the same height.

Folding and Unfolding

- To fold the walker, remove the bag, then press and fold both parts of the crossbar (Fig.5). Then hold both handles of the walker and press them inward.
- To unfold the assembled walker, hold both handles and pull them outward. Then unfold both parts of the crossbar until it locks into place.

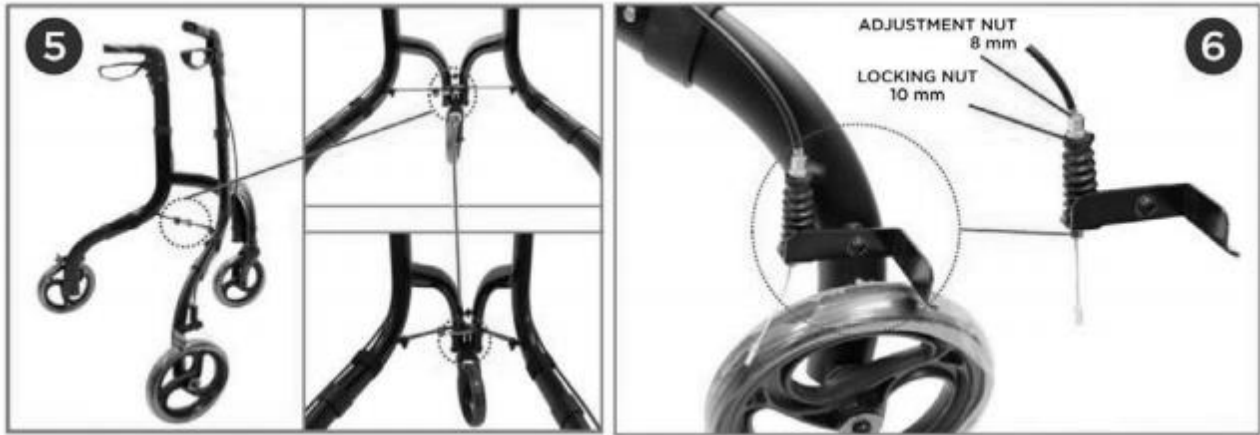

Brakes Adjustment

Use a 10mm spanner to loosen the locking nut. Then use an 8mm spanner to adjust the brake tension (there must be a small gap between the brake mechanism and the wheel, when the brake lever is released). Tighten the 10mm locking nut - make sure the adjustment nut stays in place while tightening.

Technical data

Total width	68.5 cm
Total length	67 cm
Total height	84-94 cm
Wheels' diameter	8"
Weight	4.4 kg
Max. capacity	136 kg
Max. bag capacity	5 kg

Notes

MobiQuip Ltd does not bear responsibility of improper use of the rollator, failure to comply safety regulations and misuse.

Maintenance

Clean the frame with damp cloth and dry thoroughly. Do not use strong polishing agents due to the damage of protective coatings. Clean the dirt from moving parts of the walker regularly (hinges, clamping screws, locks).

Storage

The walker should not be exposed to direct sunlight, water or moisture.

Labels

	Reference number		Manufacturer
	LOT number		Manufacture date
	Serial number		Medical device
	Note		Please read the instruction
	The manufacturer checked for compliance with essential medical devices requirements.		



MobiQuip Ltd
The Tower, Pera Business Park,
Nottingham Rd,
Melton Mowbray
LE13 0PB



Zhenjiang U-Nurse Medical Equipment
co, Ltd, Building 2, Central R&D Zone,
Science and Technology Industrial Park
Development Zone, No. 99 Jing 15th
Road, New District, Zhenjiang, Jiangsu
Province, China



IMPORTER FOR UK: (For UK customers)
Van Os Medical UK Limited
Excel House, Ashbrooke Park, Sherburn-
In-Elmet, North Yorkshire, LS25 6PJ,
United Kingdom



IMPORTER FOR EU27: (For EU27 customers)
Van Os Medical B.V.
Koperslagerij 3,
4651 SK, Steenberghe,
The Netherlands



+ 44 (0) 1977 681 400



sales@vanosmedical.com



www.vanosmedical.co.uk



+ 31 (0) 167 57 30 20



info@vanosmedical.nl



www.vanosmedical.nl

UKRP

SUNGO UK 3rd Floor, 70
Gracechurch Street,
London, EC3V 0HR,
United Kingdom

EC REP

Lotus NL B.V, Koningin
Juliaaplein 10,
Le Verd, 2595AA, The Hague,
Netherlands Amsterdam