
# MQ XL Stroller

[www.mobiquip.com](http://www.mobiquip.com)

**MQ**  
**MobiQuip**
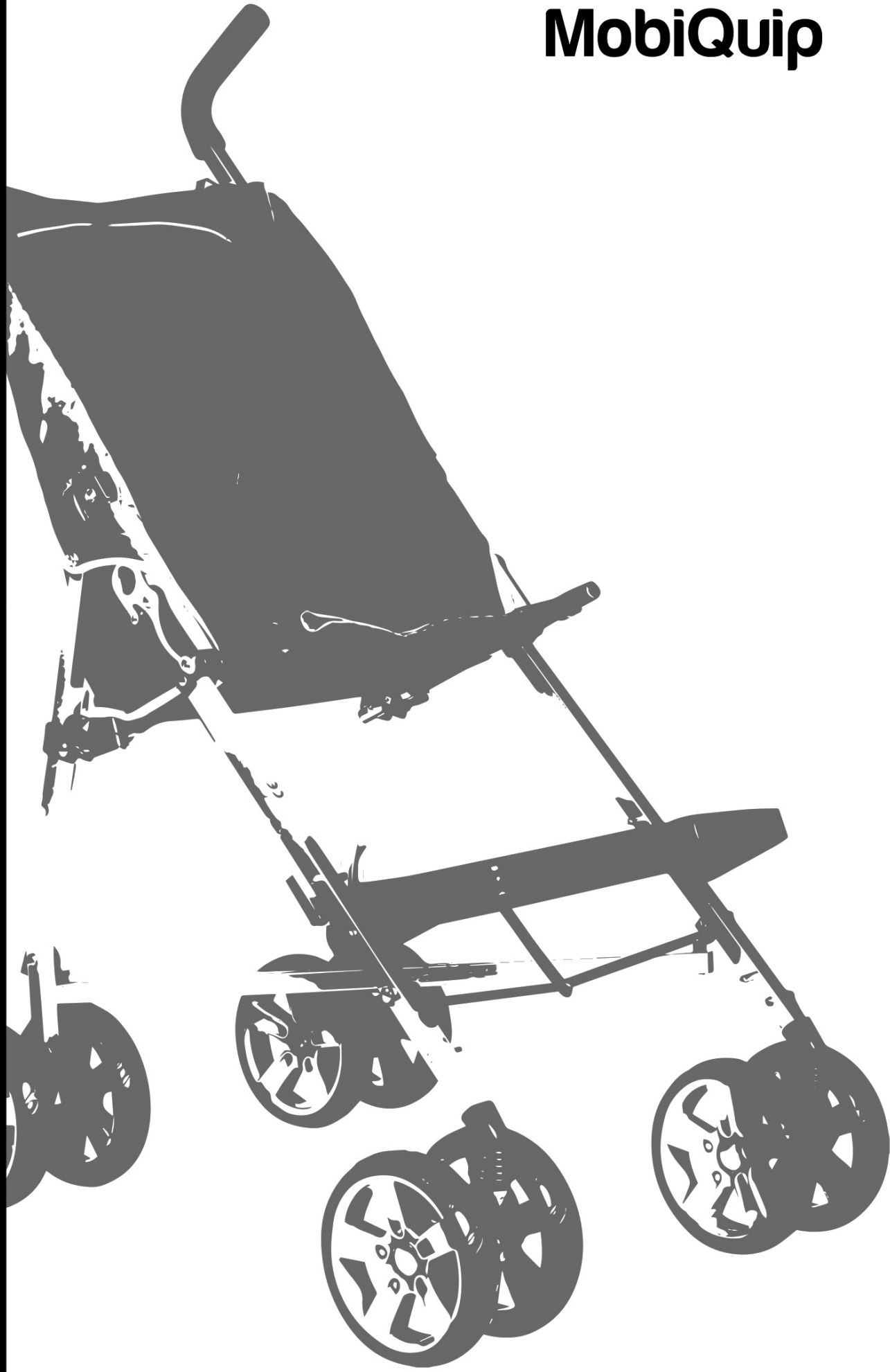
## **Table of Contents**

1. Features
2. Warning
3. Measurement
4. Assemble
5. Function

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A circular icon containing the letters 'EN' in a bold, sans-serif font.

## **WARNING!**

### **READ AND FOLLOW ALL WARNINGS AND INSTRUCTIONS! FAILURE TO DO SO COULD RESULT IN SERIOUS INJURY OR DEATH!**

NEVER leave child unattended. ALWAYS keep child in view while in the stroller.  
Avoid serious injury from falling or sliding out.  
ALWAYS use harness.  
When making adjustments to the stroller, ALWAYS ensure that all parts of the child's body are clear of any moving parts of the stroller.  
DO NOT use stroller if any parts are missing or broken.  
DO NOT allow child to climb on or around stroller.  
DO NOT allow child to stand on the footrest. Stroller may tip over and injure child.  
NEVER leave stroller on a hill or incline, even with brakes set, stroller may tip injuring the child.  
To avoid accidental rolling, ALWAYS set brakes when stroller is not being pushed.

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EN

| Measurement         |           |
|---------------------|-----------|
| Width               | 55cm      |
| Height              | 112cm     |
| Length              | 96cm      |
| Footrest height     | 112cm     |
| Footrest dimensions | 44 x 19cm |
| Seat width          | 40cm      |
| Seat depth          | 27cm      |
| Back height         | 56cm      |
| Back width          | 40cm      |
| Product weight      | 6.9kg     |
| Weight capacity     | 50kg      |

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Plegar el cochecito

**EN****Product Identification**

Your stroller is equipped with a frame label, which is located on the frame (under the seat) of your stroller. The frame label contains information about your stroller, for example the product name, manufacturer name, the name of the authorised representative(s), Unique Device Identifier (UDI) and technical data of your stroller. Below are example frame labels, the information on the frame label of your stroller will correspond to your specific stroller. Van Os Medical UK Limited is the importer of this product for the UK. This product is manufactured by Lerado Zhongshan Peaceful Cove Business Trading Co., Ltd, who are represented in the United Kingdom by UK APLAN CORPORATION LTD. Mobiquip Limited is the exclusive distributor of this mobility product. The manufacturer's product code for your stroller is MK1000.

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### Mobiquip XL Stroller (MK1000)

Paediatric travel buggy for indoor and outdoor use

03020160MQ

Lerado Zhongshan Peaceful Cove  
Business Trading Co., Ltd.  
No. 26, Guangfu Road, Dongsheng Town,  
Guangdong Province, Zhongshan City,  
528414, China

IMPORTER FOR UK:  
Van Os Medical UK Limited  
Excel House, Ashbrooke Park,  
Sherburn-In-Elmet, North Yorkshire  
LS25 6PJ

UK APLAN CORPORATION LTD  
Wandle House Associates,  
47 Wandle House, Wandle Road, Croydon, Surrey,  
CR0 1DF, United Kingdom

22-07-2021

22-07-2025

+

\* 1 2 3 4 0 1 2 3 \*

\* 0 1 0 4 8 9 7 1 1 9 8 9 0 0 2 0 1 1 2 1 0 7 2 2 1 7 2 5 0 7 2 2 1 0 1 2 3 4 2 1 0 1 2 3 \*

### Product Identification

Van Os Medical B.V. is the importer of this product for the EU27. This product is manufactured by Lerado Zhongshan Peaceful Cove Business Trading Co., Ltd, who are represented in the European Union by Y. SUNGHANDELSVERTRETUNG. Mobiquip Limited is the exclusive distributor of this mobility product. The manufacturer's product code for your stroller is MK1000.

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## Mobiquip XL Stroller (MK1000)

Paediatric travel buggy for indoor and outdoor use

03020160MQ

Lerado Zhongshan Peaceful Cove  
Business Trading Co., Ltd.  
No. 26 Guangfu Road, Dongsheng Town,  
Guangdong Province, Zhongshan City,  
528414, China

### IMPORTER FOR EU27:

Van Os Medical B.V.  
Koperslagerij 3,  
4651 SK, Steenbergse,  
The Netherlands

Y. SUNG HANDELSVERTRETUNG

Düsselthaler Str. 24,  
40211 Düsseldorf,  
Germany

22-07-2021

22-07-2025

+

\* 1 2 3 4 0 1 2 3 \*

\* 0 1 0 4 8 9 7 1 1 9 8 9 0 0 2 0 1 1 2 1 0 7 2 2 1 7 2 5 0 7 2 2 1 0 1 2 3 4 2 1 0 1 2 3 \*

MobiQuip Ltd  
The Tower, Pera Business Park,  
Nottingham Rd,  
Melton Mowbray  
LE13 0PB

Excel Mobility Manufacturing Co. Ltd.  
Rm 1907, 19/F, Dominion Centre,  
43-59 Queen's Road East, Wan Chai  
Hong Kong

IMPORTER FOR UK: (For UK customers)  
Van Os Medical UK Limited  
Excel House, Ashbrooke Park, Sherburn-  
In-Elmet, North Yorkshire, LS25 6PJ,  
United Kingdom

IMPORTER FOR EU27: (For EU27 customers)  
Van Os Medical B.V.  
Koperslagerij 3,  
4651 SK, Steenberghe,  
The Netherlands

+ 44 (0) 1977 681 400

[sales@vanosmedical.com](mailto:sales@vanosmedical.com)

[www.vanosmedical.co.uk](http://www.vanosmedical.co.uk)

+ 31 (0) 167 57 30 20

[info@vanosmedical.nl](mailto:info@vanosmedical.nl)

[www.vanosmedical.nl](http://www.vanosmedical.nl)

**UKRP**

SUNGO UK  
3rd Floor, 70 Gracechurch Street,  
London, EC3V 0HR,  
United Kingdom

**EC REP**

SUNGO EUROPE B.V.  
Olympisch Stadion 24, 1076 DE,  
The Netherlands